
nedis

Stomach and Back Heating Pad

with 6 heat levels

PEHP112CGY

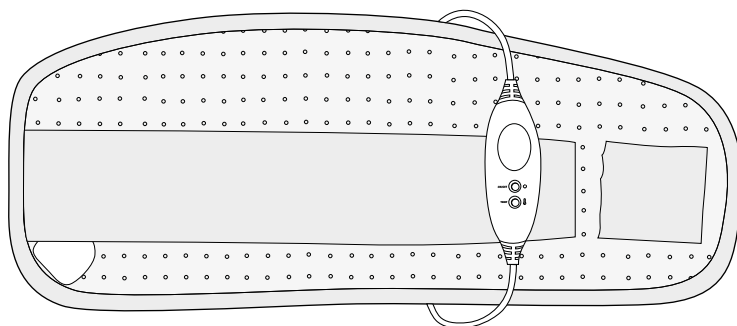
ned.is/pehp112cgy

| | | |
|----|---------------------------|----|
| EN | Quick start guide | 3 |
| DE | Kurzanleitung | 4 |
| FR | Guide de démarrage rapide | 6 |
| NL | Snelstartgids | 8 |
| IT | Guida rapida all'avvio | 9 |
| ES | Guía de inicio rápido | 11 |
| PT | Guia de iniciação rápida | 12 |
| SV | Snabbstartsguide | 14 |
| FI | Pika-aloitusopas | 16 |

| | | |
|----|-----------------------------|----|
| NO | Hurtigguide | 17 |
| DA | Vejledning til hurtig start | 19 |
| HU | Gyors beüzemelési útmutató | 20 |
| PL | Przewodnik Szybki start | 22 |
| EL | Οδηγός γρήγορης εκκίνησης | 23 |
| SK | Rýchly návod | 25 |
| CS | Rychlý návod | 27 |
| RO | Ghid rapid de inițiere | 28 |

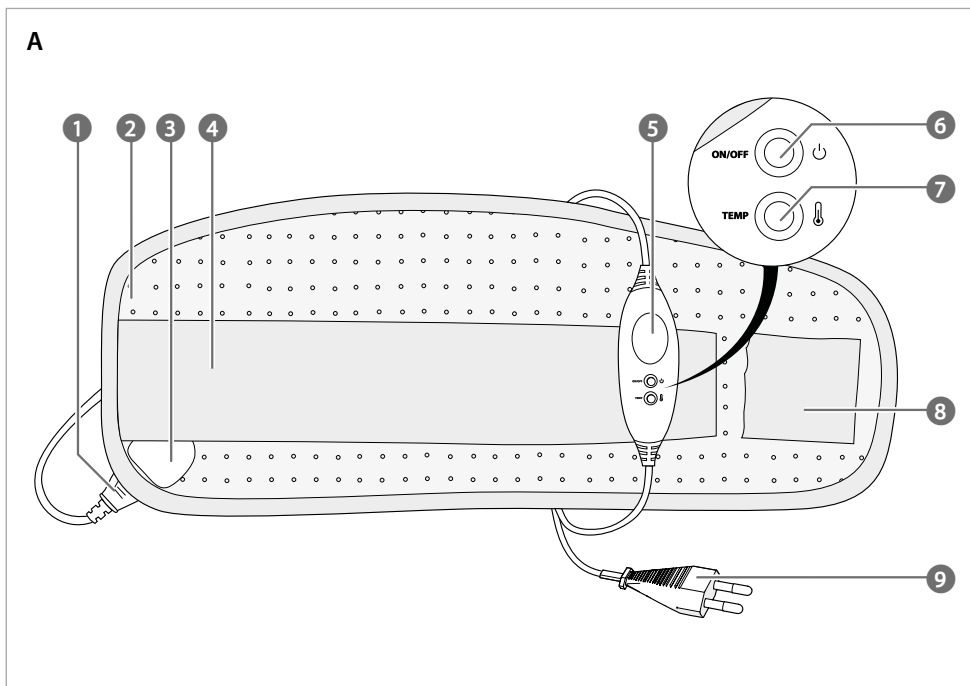Specifications

| | |
|--------------------------|------------------------------|
| Product | Stomach and Back Heating Pad |
| Article number | PEHP112CGY |
| Dimensions (l x w x h) | 690 x 280 mm |
| Power supply | 220 - 240 V ~ 50 Hz |
| Power consumption | 100 W |
| Max. temperature setting | Up to 50 °C |
| Material | Polyester and sponge |

Stomach and Back Heating Pad

PEHP112CGY

For more information see the extended manual online:
ned.is/pehp112cgy

Intended use

The Nedis PEHP112CGY is a stomach and back heating pad with a detachable controller that includes a timer, 6 heat levels and LED display.

The product offers protection against overheating and is machine washable.

The product is intended for indoor use only.

This product can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the product in a safe way and understand the hazards involved. Children shall not play with the product. Cleaning and user maintenance shall not be done by children without supervision.

This product is not intended for persons insensitive to heat and other vulnerable persons, such as children under the age of 3, who are unable to react to overheating.

The product is not intended for professional use.

The product is not intended for medical use in hospitals.

This product is not intended for warming animals.

Any modification of the product may have consequences for safety, warranty and proper functioning.

Main parts (image A)

| | | | |
|----|--------------|---|---------------------|
| ! | Power plug | & | On/off button |
| " | Heating Pad | ' | Temperature button |
| # | Connector | (| Hook and loop strip |
| \$ | Elastic band |) | Power cable |
| % | LED display | | |

Safety instructions

⚠ WARNING

- Only use the product as described in this document.
- Ensure you have fully read and understood the instructions in this document before you install or use the product. Keep the packaging and this document for future reference.
- Do not use the product if a part is damaged or defective. Replace a damaged or defective product immediately.
- This product may only be serviced by a qualified technician for maintenance to reduce the risk of electric shock.
- Unplug the product from the power source if problems occur.
- Do not unplug the product by pulling on the cable. Always grasp the plug and pull.
- Make sure the product and power cable do not get trapped when in use.
- Do not let the power cable touch hot surfaces.
- Unplug the product when not in use and before cleaning.
- Do not drop the product and avoid bumping.
- Always use the voltage (A.C. Only) marked on the rating label.
- Only use the product with the provided controller.
- Always switch off the controller before you insert or remove a plug.
- Do not expose the product to water, rain, moisture or high humidity.
- Do not immerse the product in water for cleaning.
- Do not use the product when it is damp or wet.
- Do not use the product in wet or damp environments.
- Do not use together with another heating appliance.
- Do not expose the product to direct sunlight, naked flames or

heat.

- Do not use the heating pad directly on your skin.
- Prolonged use at the highest temperature level could lead to burns.
- It is recommended to adjust the temperature to a setting for continuous use when the product is likely to be used for an extended period.
- Do not fall asleep when using the product.
- Children should be supervised to ensure that they do not play with the product.
- Cleaning and user maintenance shall not be done by children without supervision.
- The cable and control unit of the product can rise the risks of entanglement, strangulation, tripping or treading if not correctly arranged. Make sure that excess ties and electric cables are arranged in a safe way.

Installing the product

1. Connect the power plug **A!** to the connector **A#**.
2. Connect the power cable **A)** to a power outlet.

Preheat

⚠ Do not use the product when it is damp or wet.

⚠ Make sure the product and power cable do not get trapped when in use.

i Preheat 10 minutes before using.

1. Press and hold the on/off button **A&** for 3 seconds.

The LED display **A%** shows temperature level 3.

2. Press the temperature button **A'** until **A%** shows temperature level 6.

Allow the product to preheat for 10 minutes. The product can now be used.

Using the product

⚠ Do not use the heating pad directly on your skin.

⚠ Prolonged use at the highest temperature level could lead to burns.

⚠ It is recommended to adjust the temperature to a setting for continuous use when the product is likely to be used for an extended period.

i Use the back of your hand to test the temperature and to avoid burn injuries to the skin.

1. Press **A'** to set the desired temperature between 1 (low) and 6 (high).

2. Close the elastic band **A\$** with **A(**.

The product automatically switches off after 90 minutes.

After use

1. Press and hold **A&** for 3 seconds to switch off the product.

2. Unplug the product from the power outlet and allow it to cool down.

3. Disconnect **A!** from **A#**.

Storing the product

1. Make sure the product is dry.

⚠ Allow the product to dry naturally. Do not dry the product in the sun, do not use a product such as a hair dryer to dry it and do not connect it to the power outlet to dry it.

2. Gently fold the product and place it back into the packaging.

i Do not crease the product.

3. Store in a cool dry place.

i Do not place any objects on the product.

Troubleshooting

| Problem | Solution |
|---|---|
| A% shows 'P'. The product probably has an internal connection problem. | <ol style="list-style-type: none"> Press and hold A& for 3 seconds to switch off the product. Check if the product is properly connected to the controller and power outlet. Press and hold A& for 3 seconds to switch on the product. If steps 1 - 3 did not give a result: <ol style="list-style-type: none"> Unplug the product from the power source. Contact a qualified technician for service/repair. |
| A% shows 'E'. The product is probably defective. | <ol style="list-style-type: none"> Unplug the product from the power source. Contact a qualified technician for service/repair. |

Cleaning

- Unplug the product from the power outlet and allow it to cool down before cleaning.
- Clean the controller regularly with a soft, clean, dry cloth. Avoid abrasives that can damage the surface.
- To clean the product, hand wiping it down with a soft damp cloth with very mild detergent, mildly process at a maximum washing temperature of 40°C. Allow it to dry thoroughly before using or storing it. When washed, drape it over a clothes line to dry.
- Do not use aggressive chemical cleaning agents such as ammonia, acid or acetone when cleaning the product.
- Do not use clothes pegs to secure it.
- Allow the product to dry naturally. Do not dry the product in the sun, do not use a product such as a hair dryer to dry it and do not connect it to the power outlet to dry it.
- Do not immerse the controller in water or get it wet.
- Do not wring out the product.
- Do not use mothproofing agents.

Cleaning in a washing machine

A! Only the heating pad **A!** is suited for machine washing. However, it is advised to handwash the product for a longer product life.

- Disconnect **A!** from **A#**.
- Place **A!** in a net bag and seal it.
- Use a small amount of detergent.
- Set the washing machine to a "Delicate Wash" with a maximum temperature of 40°C.
- When washed, remove **A!** from the net bag.
- Drape the product over a clothes line to dry.

A! Do not use clothes pegs to secure it.

| | | | |
|--|--|---|-----------------------------------|
|  | Read the instructions before use. |  | Do not tumble dry. |
|  | Do not insert pins. |  | Do not iron. |
|  | Do not use folded or rucked. |  | Do not dry clean. |
|  | Not to be used by very young children (0-3 years). |  | Maximum washing temperature 40°C. |
|  | Do not bleach. | | |

Maintenance

Unplug the product from the power outlet and allow it to cool down before maintenance.

Do not try to open or repair the product, controller or cable. It must be serviced by a qualified technician.

Examine the product frequently to check that there are no signs of wear or damage. If any damage/wear is spotted, do not use the product and inquire advice from a product specialist.

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